After-sales service

Thank you for purchasing everbright's water flosser, please note the following matters:

1. Request the seller to issue effective sale evidence to be the proof of maintenance during your purchase.
2. The product is freely maintained and repaired due to non-artificial reason within one year in accordance with related regulations of the State (based on the invoice date of purchasing);
3. If the following situation occurs within the warranty period, certain maintenance cost will be charged:
 - ①. Damage caused by consumer's misuse, incorrect maintenance and improper management;
 - ②. Damage caused by any unauthorized repairman;
 - ③. There is no relevant receipt or invoice, and the product has exceeded the production date for one year;
 - ④. Damage caused by using for business purpose.
 - ⑤. Damage caused by force majeure (such as earthquake and fire).
4. Please carefully complete the following table, and send it to the manufacturer's designated maintenance site, and attach the product with the repair requirement

Maintenance Records

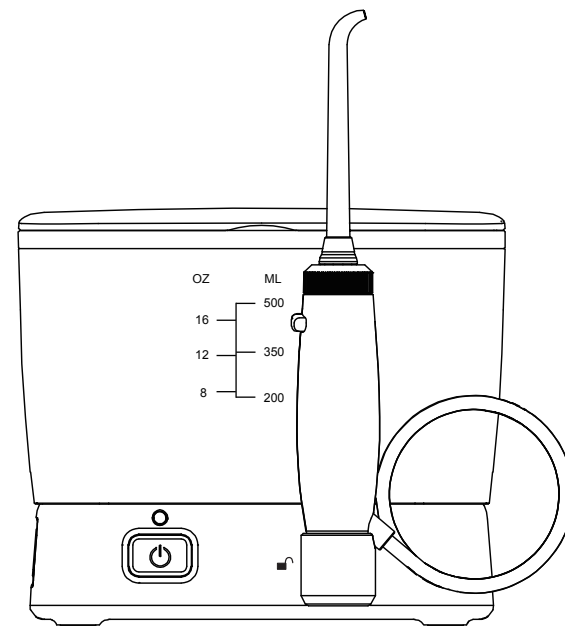
Name of customer		Tel.	
Product type		Product color	
Date of purchase		Where to buy? (seller)	
Date of maintenance			
Delivery address			

Tandenborstel.com B.V.
Heusdensebaan 50, 5061 PS Oisterwijk, Netherlands
Tel: 31 (0) 88 – 52 51 500
Web: www.tandenborstel.com

ever*bright.*

Water Flosser

Instruction Manual

Model: FC2580

UL standard IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not immerse in water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
- Do not reach for a product that has fallen into water.
Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.

To reduce the risk of burns, electric shock, fire, serious injury or damage issues:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not use this product if it has a broken cord or plug. In case of not working properly, or it has been dropped into any liquid, contact your local distributor.
- Do not direct water flow under the tongue, into the ear, nose or other delicate organs. This product is capable of producing pressures that may cause serious damage to these organs. Refer to operating instructions for correct usage.
- Use this product only as indicated in these instructions or suggestions by your dental professional. It's normal phenomenon if gum bleed a little when using water flosser for the first time. If bleed much, please stop using and ask for your dentist advice.
- Fill reservoir with water or other solutions recommended by dental professional.
- Do not drop or insert any foreign object into the product.
Keep charger and cord away from heated surfaces.
- Do not operate where oxygen or aerosol sprays are being used.
- Do not use mouth washes/rinses that contain iodine, bleach, or tea tree oil.
- Do not use while wearing any oral jewelry. Please put off prior to use.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advise you to receive antibiotic premedication, you should consult your dentist before using this instrument or any other oral hygiene aid.
- For household use only. Do not use outdoors.
- Closely instruct and supervise children and individuals with special needs in the proper use of this product.
- Children should be supervised to ensure that they do not play with the appliance.

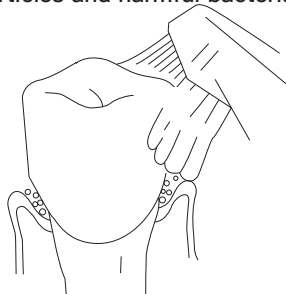
Thank you for purchasing our Water Flosser!

Welcome to use "everbright" Water Flosser, it would bring you healthier teeth and brighter smile!


The water flosser use pulse flow to clean teeth, which could clean harmful bacteria that deeply hidden under gum, inter-dental -- these area cannot be cleaned by only brushing your teeth.

Traditional oral care:

Traditional toothbrush cannot clean away deeply hidden food particles and harmful bacteria.

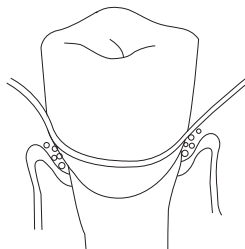
The toothbrush could only touch limited location under the gum line.

Everbright Water Water Flosser

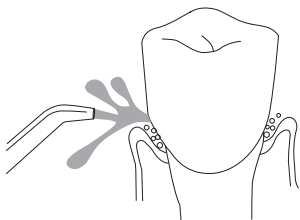It could take away the harmful bacteria that being deeply hidden under and the gum lines the inter-dental.


Catalog

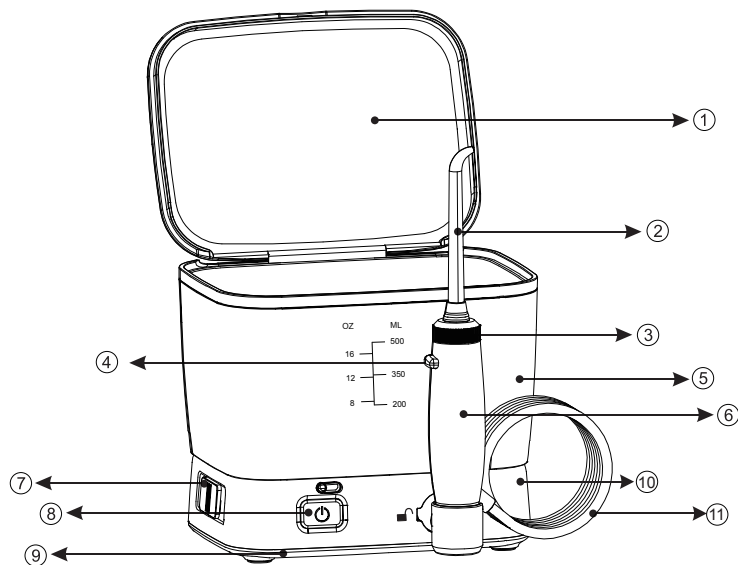- UL standard----- 1
- Product structure ----- 3
- Product accessories ----- 4
- Instructions ----- 5
- Inserting and removing the tip ----- 6
- Adjusting water pressure ----- 6
- Using Water Flosser ----- 7
- Connecting the base and the handle----- 7
- Working time and charging----- 8
- Installation and un-installation----- 8
- Using tips----- 9
- Pause----- 10
- Common trouble shooting----- 12

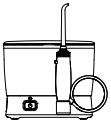
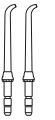
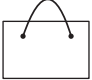
The Company reserves the right of final interpretation for any misprint or misunderstanding of the operation manual. In case of product or information update will be included in the new manual without prior notice; if there is any change in the color of the product, the actual goods in package shall prevail.

Product structure

- | | |
|--------------------------|-------------------------|
| 1. Water tank lid | 2. Jet Tip |
| 3. Tip Rotating knob | 4. Jet Tip Eject Button |
| 5. Water tank | 6. Handle |
| 7. Water pressure switch | 8. Power |
| 9. Base | 10. Body |
| 11. Water pipe | |


Product accessories

S/N	Name	Quantity	Picture
1	Main Device	1	
2	Standard Jet tip	2	
3	Charging cable	1	
4	Instruction manual	1	
5	Handle	1	
6	Bag	1	

Because gingiva condition differs from people to people, the first use of the Water Flosser may have gum bleeding symptom. This is normal. Please feel free to use it! However, if gum bleeding too much, and last long time, please advise your dentist.


Instructions

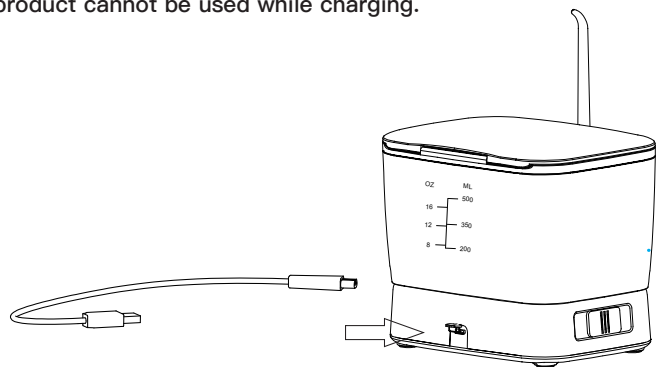IMPORTANT! Please fully charge the Water Flosser for the first use by plugging the USB cable into the charging socket with power. The product cannot be used while charging.

Filling water

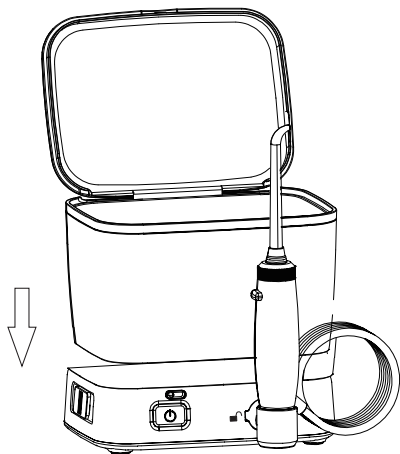
Remove the water tank and water tank lid from the base. Fill pure water into the water tank, place it on the base and lightly press it down.

7


Inserting and removing the tip

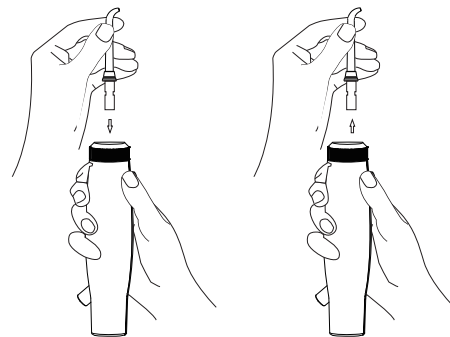Insert the jet tip in the hole of the tip knob top, make sure that the jet tip is inserted in position till a clicking sound. If needing to remove the tip from the handle, press the eject button while pull off jet tip from the handle.

Adjusting water pressure

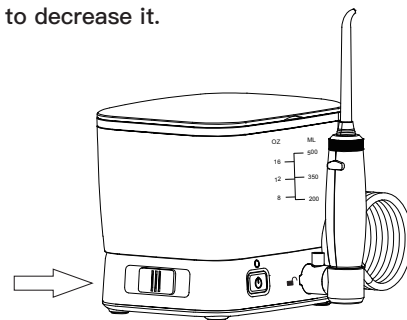Push forward the button to increase water pressure, push back to decrease it.


8

Using Water Flosser

Bend your body to the washbasin, put the jet tip into your mouth to aim it to your teeth, press the switch with the other hand and set your preferred pressure.

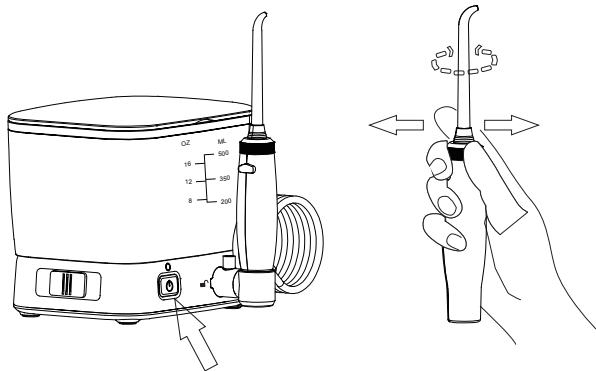
Rotating the knob at the top of the handle to change water flow direction.

Connecting the base and the handle

Insert the end of the handle to the outlet of the device base as marking, then rotate 90° in clockwise direction to erect the handle. There shall not have leakage from the connection.

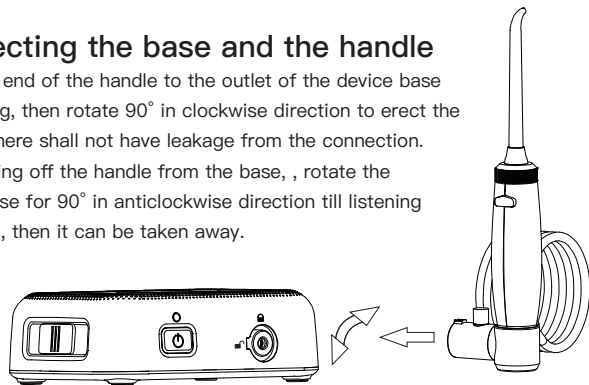When taking off the handle from the base, rotate the handle base for 90° in anticlockwise direction till listening the sound, then it can be taken away.


Working time and charging

The regular working time of the product is 150s, if you need to use it for more time, turn on power again. It could be stopped by pressing the switch during working.

When low battery, the white indicator will flash quickly. There will be a white indicator flashing slowly during charging, when the indicator is lighted for a long time, it means full charge. The product cannot be used during charging.

Installation and un-installation



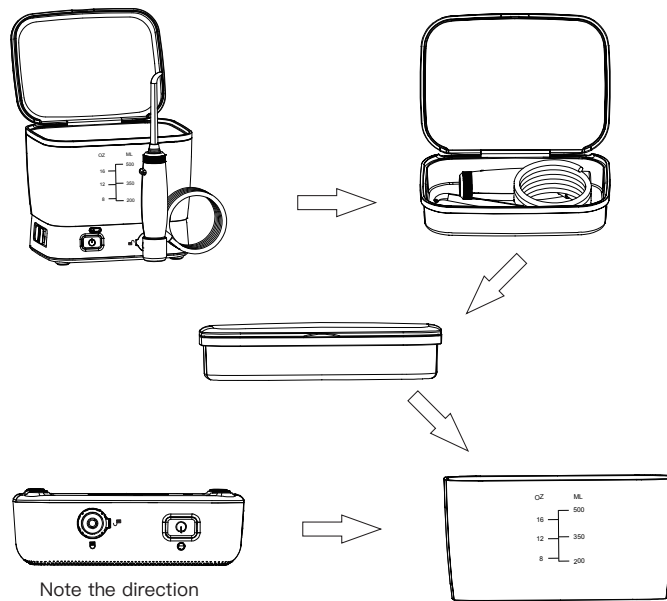
The water tank lid can accommodate the handle and jet tip, the tank can contain device base, then cover water tank lid.

Note the direction


Using tips

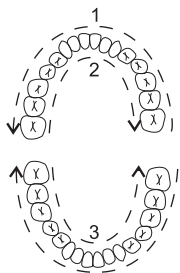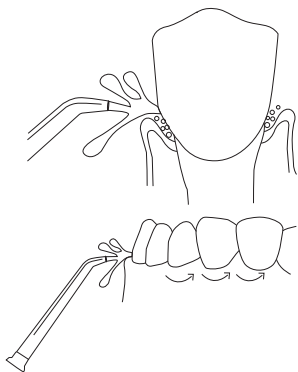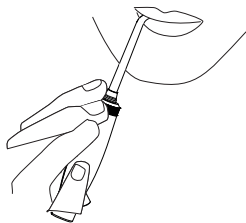Aiming at the gingiva by the water flow in 90° direction, slightly close your mouth to prevent the water splash everywhere, but allow the water to flow along the corner of your mouth to the basin. For better cleaning, you should start the cleaning from your molar tooth to incisor. Lightly move the jet tip along the gingiva and stay for a while. Repeat the same operation till everywhere of your teeth are cleaned.

Pause

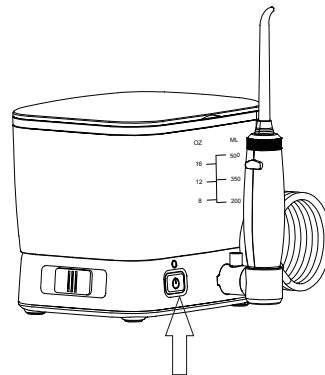
You could stop the machine by pressing the power button or turning the pressure pushing to minimum value.

Instruction of the tip

Standard Jet tip

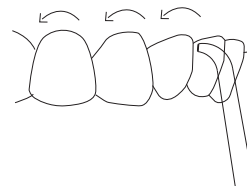360° clean. And the place which cannot be cleaned by the toothbrush. It is more effective than the dental floss and is healthier for your gingiva.

Cleaning:

If any specific solution is stored in the water tank, take little time to clean the water tank after using.

Fill some warm(not hot) water into the water tank and point jet tip to the washbasin. Switch on the machine to run it till the water is emptied. This helps to prevent the blocking and increase the life of the product.

When there is necessary for washing, clean the product with soft cloth and neutral wear-free detergent.

Maintenance

- 1、 After using the Water Flosser every time, invert the water tank, till no water or cleaner and dried, and place it in cool and dry place with cleaned device. If the machine is not used for a long time, it should be put into the package and store well.
- 2、 When there is something wrong with the machine, you could consult the retailer.
- 3、 Please do not dismantle or try to disassemble the product by yourself, and the failure caused by this will not be included the warranty.
- 4、 If the charge cable is broken, please repair it at the repair shop designated by the Company or the shop with relevant repair qualification.
- 5、 We still supply service for the Water Flosser that has exceeded the warranty period, but material cost and little manual service fee will be charged.

Common trouble shooting

Error phenomena	Reason	Treatment method
Leakage between the jet tip and handle	The jet tip is not installed in place	Remove the tip and install it in right place
Leakage from the water pipe	Damaged pipe	Replace it through the seller or repair it by the customer service
Water tank leakage	It might be no valve at the bottom of tank, or the valve is mounted in an opposite way	Replace the valve or mounted in the right direction.
Insufficient pressure	Incorrect location for the installation of the water tank	Mount the water tank again with downward pressure
Failure to start	Motor failure	Ensure whether the power supply is switched on; contact with the sales repair

The Company assures to the original consumers and owners of this new product that the product has no defects in material and process, and the warranty period is one year from the date of purchase.

Thank you once again for purchasing everbright Water Flosser!

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