

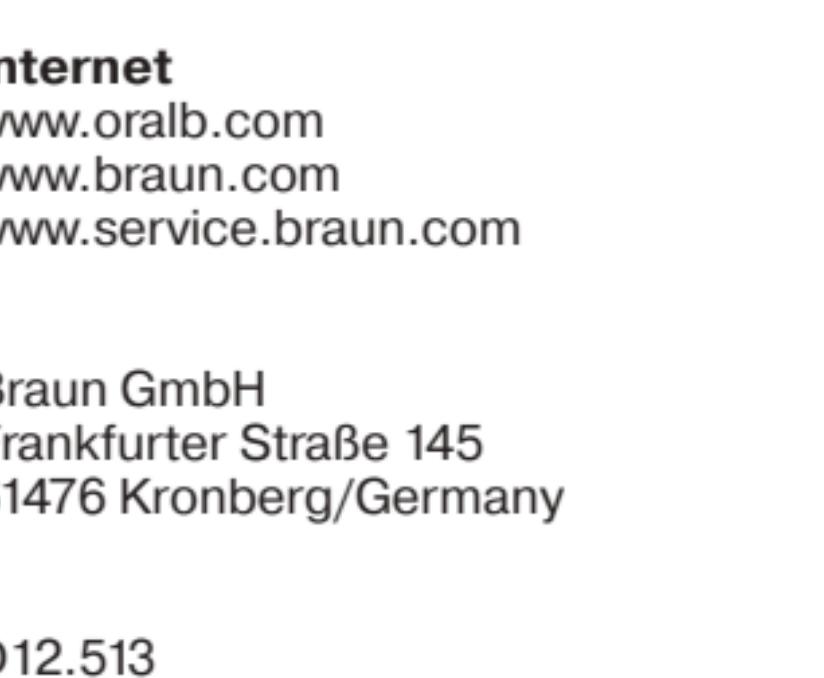
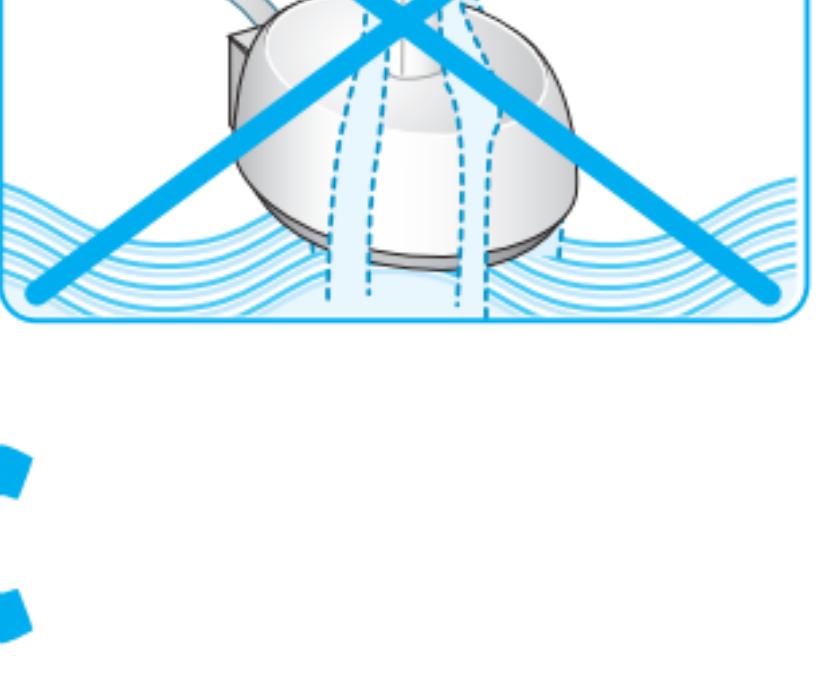
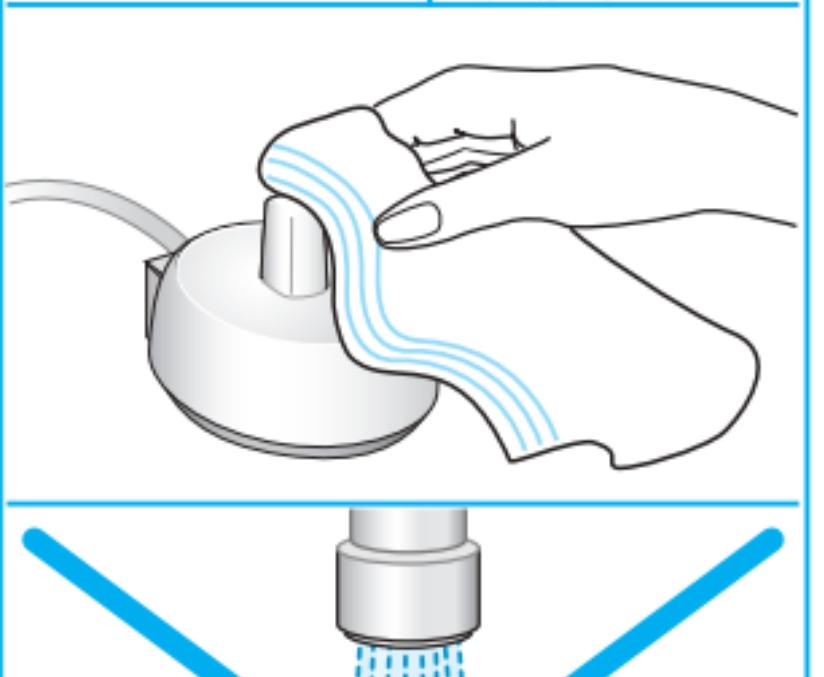
Oral-B

TRI ZONE powered by BRAUN

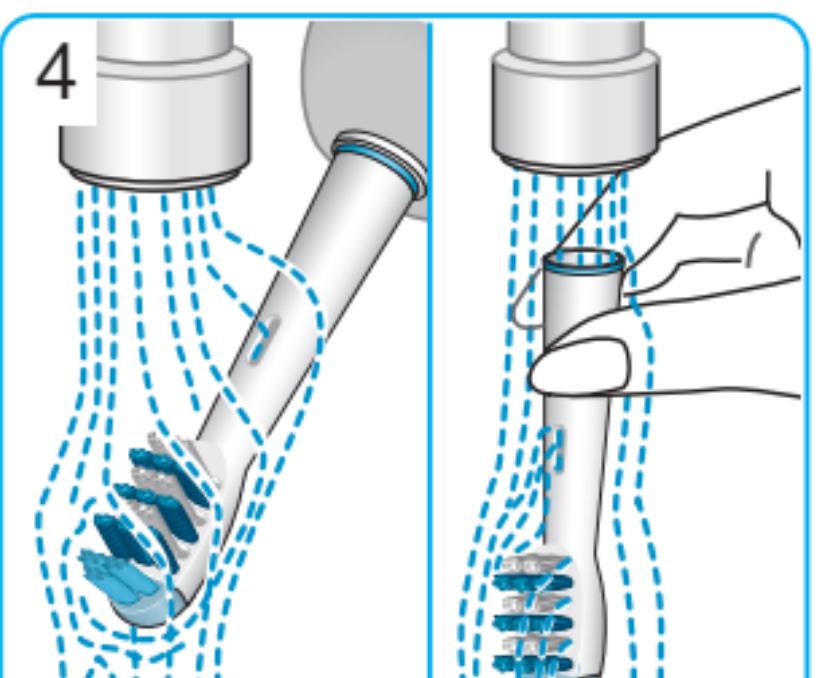Type 3757

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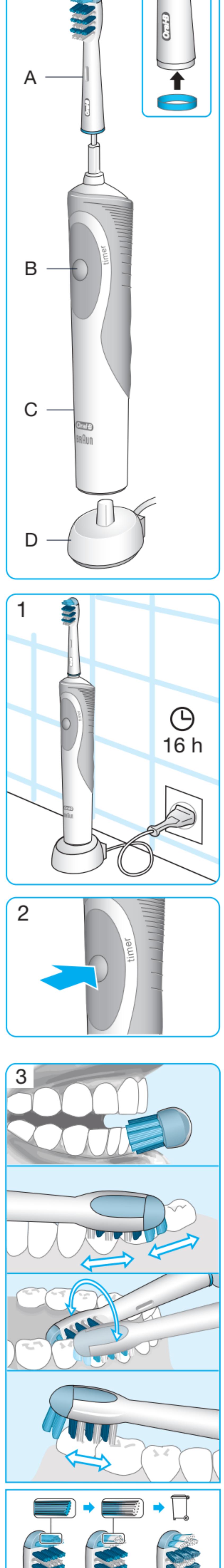plugged-in charging unit to maintain it at full power. Overcharging is impossible.

- However, for environmental reasons, Oral-B recommends to unplug the charging unit until the next recharge is required.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

Using your toothbrush

Brushing technique

This brush is designed to work with all common brushing techniques. You can continue to use your current brushing technique, e.g. you can brush just like you would do with a manual toothbrush. You also may consult your dentist or dental hygienist about the right technique for you.

We recommend the following:

Wet brush head and apply toothpaste (before switching the appliance on); you can use any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on. You start and stop brushing by pressing the on/off button of the handle (picture 2).

Place the toothbrush bristles against the teeth at a slight angle towards the gumline.

Apply light pressure and start brushing in back and forth motions. Brush first the outsides, then the insides. Finally, brush the chewing surfaces with the same back and forth motion throughout the mouth (picture 3). Brush for 2 minutes. Please ensure that you brush all four quadrants of your mouth equally.


Timer (depending on model)

An incorporated timer memorises the elapsed brushing time. A short stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

Brush heads


Oral-B TriZone brush head

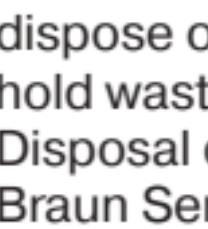
features triple zone cleaning action for outstanding plaque removal, even in between teeth.

Oral-B also offers you a variety of different Oral-B brush heads that fit your Oral-B toothbrush handle and can be used for tooth-by-tooth precision cleaning.


Oral-B FlossAction brush head

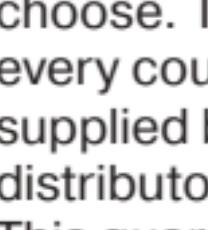
features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

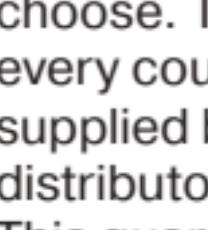cups each tooth with a curved bristle trim and reaches deep between teeth at the same time.

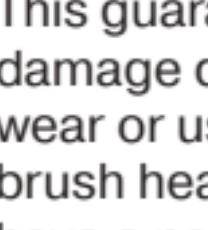features an extra soft bristle structure that is gentle on teeth and gums.

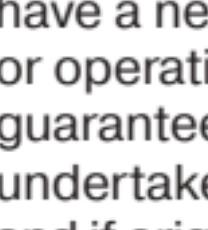features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B 3DWhite brush head.

Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums (picture 3).

We do not recommend using the Oral-B FlossAction or Oral-B 3DWhite brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

Cleaning recommendations

After brushing, rinse your brush head for several seconds under running water with the handle switched on.


Switch off the handle and remove the brush head. Rinse both parts separately under running water then wipe them dry before reassembling the toothbrush. From time to time also clean the charging unit with a damp cloth; never place the charging unit in water (picture 4).

Subject to change without notice.

This product contains rechargeable batteries.

In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life.

Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

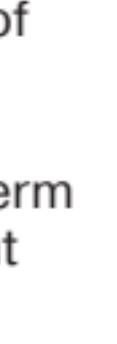 This guarantee does not cover:

damage due to improper use, normal wear or use, especially regarding the brush heads, as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

Please also read our additional warranty statements: **Warranty replacement brush heads** on page 5.

To obtain service within the guarantee period, hand in or send the complete

appliance with your sales receipt to an authorised Oral-B Braun Service Centre.

This guarantee in no way affects your rights under statutory law.

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